Surgical System

Surgic Pro2

OPERATION
MANUAL

CE 0197 MADE IN JAPAN

English	1
Deutsch.....	61
Français.....	125
Español	189
Italiano.....	251
Português.....	315

TABLE OF CONTENTS

English

1	Preface	3
1.1	User and intended purpose	3
1.2	Related manuals	3
1.3	Precautions for handling and operation	3
1.4	Symbol.....	6
1.5	Disposing of products	8
1.6	Warranty.....	8
2	Product Description	9
2.1	Package contents	9
2.2	Applicable products	10
2.3	Part names.....	11
2.4	Control panel	13
3	Preparation for Use	16
3.1	Connecting AC power cord	16
3.2	Installing foot control batteries	16
3.3	Installing foot control hanger.....	17
3.4	Mounting coolant solution hanger post	17
3.5	Installing irrigation tube.....	18
3.6	Inserting irrigation tube onto and from bag / bottle.....	19
3.7	Connecting motor cord	19
3.8	Connecting handpiece.....	20
3.9	Irrigation nozzle attachment.....	20
3.10	Attaching tube holder	21
3.11	Check before treatment.....	21
3.11.1	Compatibility check of internal spray nozzle/drill	22
4	Operation Procedure	23
4.1	Calibration function	23
4.1.1	Calibration preparation	23
4.1.2	Conducting calibration.....	24
4.2	Starting operation.....	26
4.3	Protection circuit.....	27
4.4	Sleep mode	27
4.5	Link function.....	28
4.5.1	Osseo integration monitoring device (Osseo 100+).....	29
4.5.2	Link Module.....	30
4.5.3	External terminal.....	31
5	Post-use Maintenance	32
5.1	Preparation	32
5.2	Cleaning and disinfection	33
5.3	Sterilization	35

5.4	Storage	36
6	Maintenance.....	37
6.1	Replace O-ring	37
6.2	Periodical maintenance checks.....	37
6.3	Foot control calibration.....	38
7	Setting.....	40
7.1	Setting key.....	40
7.1.1	Pairing with foot control	43
7.1.2	Foot control button guidance.....	44
7.2	Implant system	45
7.3	Programming system operation	46
7.3.1	Program skip function	47
8	Troubleshooting	48
8.1	Error code.....	48
8.2	Problems and solutions.....	49
9	Spare Parts and Option List	53
10	Specification.....	55
10.1	System specifications	55
10.2	Bluetooth specifications.....	56
10.3	Classification of equipment	56
10.4	Operation principle	56
10.5	Factory settings (implant system).....	56
11	EMC Information (Electromagnetic Compatibility Information)	58

1 Preface

In this operation manual, Surgic Pro2 is called "this product".

1.1 User and intended purpose

- User
Qualified professionals (Dentists and oral surgeons).
- Intended purpose
The dental and oral surgical equipment is intended for the following application(s):
Cutting, perforation, removal and preparation of jawbones, placement, tightening, and removal of dental implants (fixtures, abutment), and removal of third molars

1.2 Related manuals

For proper use with the combined items, refer to the following operation manuals:

Operation manual title	Operation manual No.
Ti-Max Z-SG45L / X-SG20L / X-SG20ML / X-SG25L / nanoSG20LS / X-SG65L / X-SG93L, S-Max SG20 / SG20M OPERATION MANUAL	OM-C0595E
Ti-Max X-DSG20L/X-DSG20/X-DSG20Lh/X-DSG20h OPERATION MANUAL	OM-DC0929MA
Link Module OPERATION MANUAL	OM-E0859E
Link Stand2 OPERATION MANUAL	OM-DE0941MA
New compatible terminals have been added	OM-DE1112MA

1.3 Precautions for handling and operation

- Read these safety cautions thoroughly before use and operate the product properly.
- These indicators are to allow you to use the product safely, prevent danger and harm to you and others. These are classified by degree of danger, damage, and seriousness. All indicators concern safety, be sure to follow them.

Class	Degree of risk
⚠ WARNING	Hazard that could result in serious injury or damage to the device if the safety instructions are not correctly followed.
⚠ CAUTION	Hazard that could result in light or moderate injury or damage to the device if the safety instructions are not correctly followed.
NOTICE	General product specification information highlighted to avoid product malfunction and performance reduction.

⚠ WARNING

- TO PREVENT ELECTRIC SHOCK, do not unplug the AC power cord with wet hands.
- TO PREVENT ELECTRIC SHOCK, be sure to prevent water or any fluid from getting on the control unit. If it gets on the control unit, wipe it off immediately.
- TO PREVENT ELECTRIC SHOCK, use an electrical outlet that is grounded.

WARNING

- The system may malfunction when used in the presence of an electromagnetic interference wave. Do not install the system in the vicinity of a device which emits electromagnetic waves. Turn off the main power switch of the control unit if a high frequency oscillation device such as an electrode knife is in the vicinity and used.
- Do not operate close to patients with cardiac pacemakers as there is a danger that it may affect the pacemaker.
- Do not use in rooms with a risk of explosion or near flammable substances. Also, do not use on or near patients under flammable anesthesia (nitrous oxide).
- Be sure to calibrate this product every time before use. Without calibration, there is a risk that the set torque and real applied torque output may differ.

CAUTION

- NSK surgical (oral surgery) handpieces conforming to ISO 3964 (EN ISO 3964) are available for use in this product.
- The accuracy of the torque monitor depends on the high precision performance of the handpiece attached to the micromotor. Actual torque value may not be displayed correctly. To guarantee handpiece torque accuracy that matches the monitor display system, USE ONLY the NSK contra angle handpieces listed in “2.2 Applicable products”. If other handpiece brands are connected, the torque output maybe incorrect.
- When operating the product always consider the safety of the patient.
- Read this Operation Manual before use to fully understand the product functions and file for future reference.
- This product needs special precautions regarding EMC and needs to be installed and put into service according to the EMC information. (Refer to “11 EMC Information (Electromagnetic Compatibility Information)”)
- Portable and mobile RF (Radio Frequency) communications equipment can affect this product. Do not use RF equipment in the vicinity of the product.
- The use of accessories, motors and cables other than those specified, with the exception of motors and cables sold by the manufacturer of this product as replacement parts for internal components, may result in increased EMISSIONS or decreased IMMUNITY of the control unit.
- If this product must be used adjacent to or stacked with other equipment, make sure that the device operates without malfunction.
- Do not use non-NSK genuine irrigation tubing.
- Do not attempt to disassemble the product nor tamper with the mechanism except as recommended by NSK in this operation manual. Doing so may result in an injury, electric shock or fire.
- Before each use, especially after a long period of disuse, inspect device operation to make sure there are no abnormalities. If you detect abnormal vibration, heat, sound, or the like before or during use, stop using immediately and request repair.
- Before use, carefully read the operation manual of the handpiece for attachment and removal of the bur/drill and irrigation methods.
- Do not drop, hit, or subject to excessive shock. Doing so could cause malfunctions.
- Do not turn over the control unit. The calibration load may drop.
- Install the cord/plug only when main power is OFF.
- Make sure each component is securely connected before turning on the power.

⚠ CAUTION

- Make sure that cords and other parts do not touch the control panel while the control unit power is on. It may result in malfunction.
- Do not bend the irrigation tube while the irrigation pump is operating. Doing so could cause the tube to break or come off.
- Do not use burs or drills like those below. It may cause the shank to break or come off.
 - Bent, deformed, unsymmetrically worn, rusty, chipped, or broken bur or drill
 - Bur or drill with damaged blade or shank
 - Bur or drill that does not meet specifications or has been modified
- Do not use burs or drills in excess of the rotation speed or torque recommended by the manufacturer.
- After each operation, be sure to lubricate and sterilize the handpiece as soon as (within 1 hour) it is cleaned. Blood coagulation can cause corrosion and rusting.
- Do not lubricate the motor. Oil could generate excessive heat and cause damage.
- The control unit can be cleaned with a moist cloth. Disconnect the power supply before cleaning.
- The control unit and the foot control cannot be sterilized by any method.
- When the system is very frequently used, consider keeping a small stock of replaceable parts.
- Do not disconnect the motor cord from the motor.
- Never use the device over 1 minute. Between use, please allow a period of 9 minutes rest time. Continuous use over a long time may lead to overheating of the handpiece or motor, which in turn may result in burns. If using the motor with the light on for more than 1 minute continuously, the motor will become hot. (The maximum surface temperature might reach 50 °C depending on by the situation in which the system is used.) In such cases, stop using the motor until the light part cools down.

Operating Mode

Non-continuous Operation	
ON	OFF
1 min	9 min

- Do not look at the LED light directly. Doing so poses a risk of eye injury.
- If any abnormalities (dark, does not light, flashes, etc.) are found with the LED, stop using it immediately and contact an authorized NSK dealer.
- For safety reasons, install the control unit in a place where the AC power cord can be easily removed. (It is possible to disconnect the control unit from the power source by removing the AC power cord.)
- Following products are delivered in a non-sterile condition and must be autoclaved prior to use. Handpiece, motor with motor cord, handpiece stand, tube holder, protection plug, calibration bur.
- Do not use parts other than the components of this product. For handpieces, use NSK surgical handpieces conforming to ISO 3964 (EN ISO 3964).
- Please report any serious incident that may occur in connection with this device and its accessories to the manufacturer as well as to your national authority.

NOTICE

- Turn off the main power switch after each use.
 - For service requirements and spare parts, contact an authorized NSK dealer.
 - No special training is required for Qualified Professionals to operate this device.
-

1.4 Symbol

	Battery cover close
	Battery cover open
	The instrument can be used with a washer-disinfector for thermal disinfection.
	This product can be sterilized in a steam sterilizer at 135°C.
	Consult operation instructions
	Caution
	Type B applied part.
	Marking on the outside of equipment or equipment parts that include RF transmitters or that apply RF electromagnetic energy for diagnosis or treatment.
 IPX8	Protected against the effects of continuous immersion in water.
	This product is designed not to become the ignition source in air and flammable anesthetic gas.
	Dispose of this device and its accessories via methods approved for electronic device and in compliance with the Directive 2012/19/EU
	Conforms to European “Medical Device Regulation”
	Mandatory conformance mark on many products placed on the market in the European Economic Area
	TUV Rheinland of North America is a Nationally Recognized Testing Laboratory (NRTL) in the United States and is accredited by the Standards Council of Canada to certify electro-medical products with Canadian National Standards.

Regulatory compliance mark (RCM) indicates compliance with electrical safety, EMC, EME and telecommunications requirements for Australia and New Zealand.

Korean KC compliance marking for wireless product

Brazil ANATEL compliance marking for wireless product

National Communications Commission (NCC)
Taiwan NCC compliance marking for wireless equipment.

Manufacturer

Importer

Authorized representative in the European Community

Authorised representative for Switzerland

Medical Device

Caution: U.S. Federal law restricts this device to sale by or on the order of a licensed physician.

Catalog number (Order code)

Serial No.

Batch code

Date of manufacture

Unique Device Identifier

GS1 DataMatrix for Unique Device Identifier.

Contains hazardous substances

Temperature limitation

Humidity limitation

Atmospheric pressure limitation

Keep dry

Fragile / handle with care

This is the correct upright position of the distribution packages for transport and/or storage.

1.5 Disposing of products

When your medical device has reached the end of its service life, it must be disposed of via a specialist collection, removal, recycling or destruction channel. Prior to do so, in order to avoid the health risks to operators handling the disposal of medical equipment, as well as the risks of environmental contamination caused thereof, you need to confirm the equipment is decontaminated.

Therefore contact your nearest NSK dealer, or one of the NSK branches to find out how to proceed. Refer to the link below:

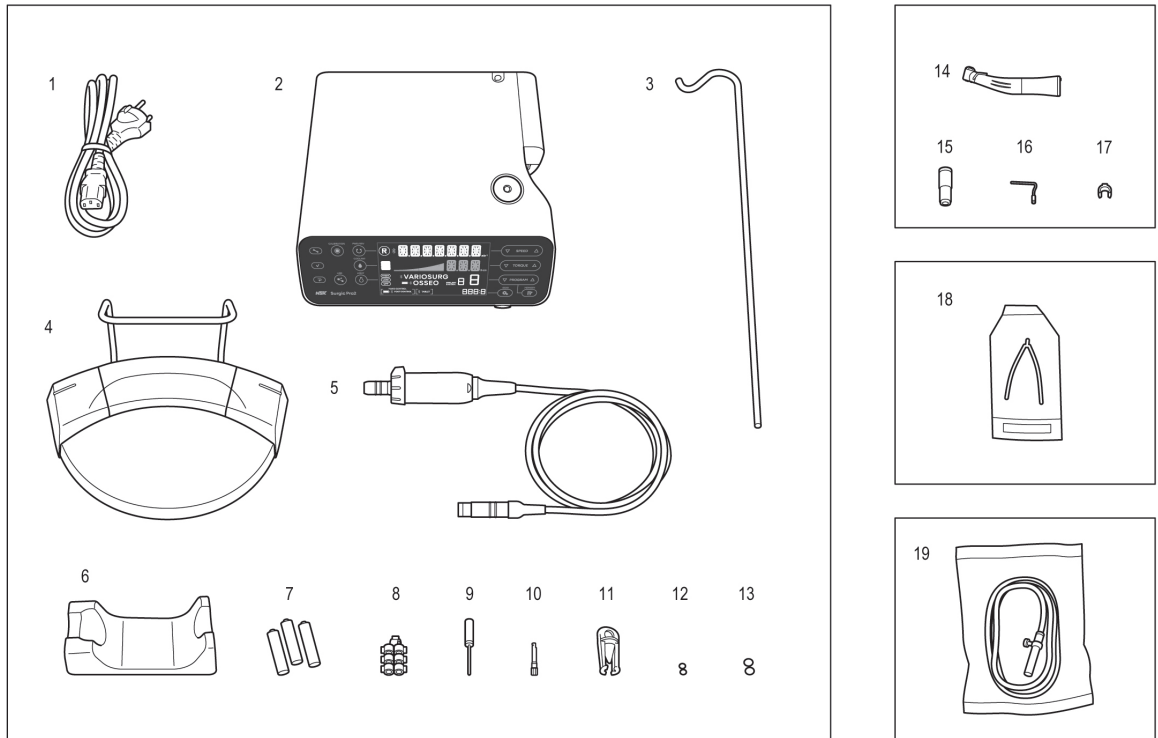
1.6 Warranty

NSK products are warranted against manufacturing errors and defects in materials. NSK reserves the right to analyze and determine the cause of any problem. Warranty is voided should the product be not used correctly or for the intended purpose or has been tampered with by unqualified personnel or has had non NSK parts installed. Replacement parts are available for seven years beyond discontinuation of the model. Contact an authorized NSK dealer if repairs are necessary.

2 Product Description

2.1 Package contents

No.	Part name	Quantity
1	AC Power Cord	1
2	Control Unit	1
3	Coolant Solution Hanger Post	1
4	Foot Control	1
5	Optic Motor / Non-Optic Motor (with motor cord) **	1
6	Handpiece Stand	1
7	AAA Alkaline Battery 3 Pack	1
8	Tube Holder	7
9	Cleaning Wire	1

** Either handpiece / motor is included.

No.	Part name	Quantity
10	Calibration Bur	1
11	Protection Plug	1
12	O-Ring (Protection Plug)	2
13	O-Ring (Motor)	2
14	Optic Handpiece / Non-Optic Handpiece **	1
15	E-Type Spray Nozzle	1
16	Internal Spray Nozzle	1
17	Nozzle Holder	1
18	Y-Connector	1
19	Irrigation Tube	3

2.2 Applicable products

Products that can be used in combination with the device

<Contra angle handpiece>

Listed below are the contra angle handpieces that should be used with the micromotor of this product to guarantee torque accuracy.

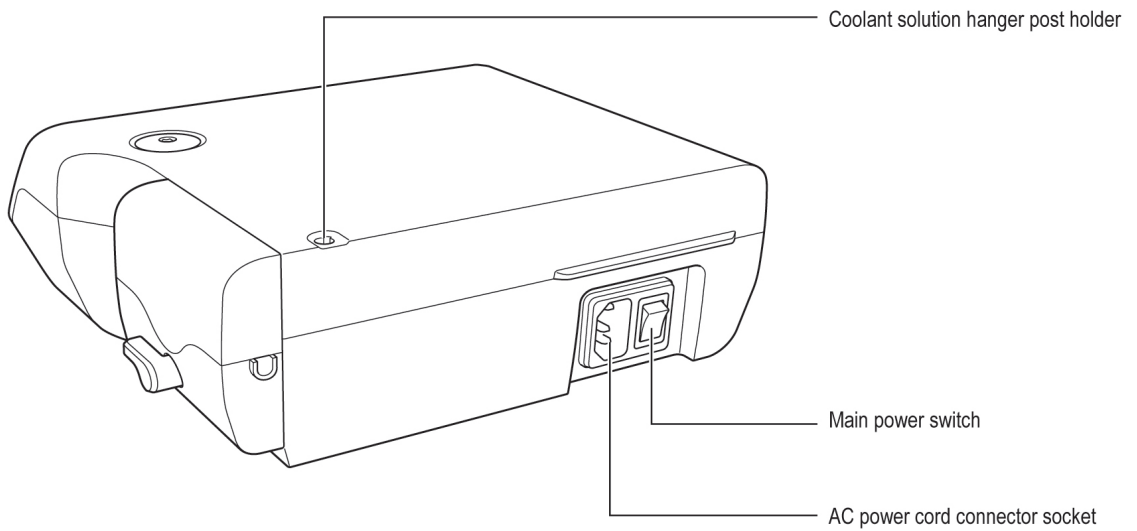
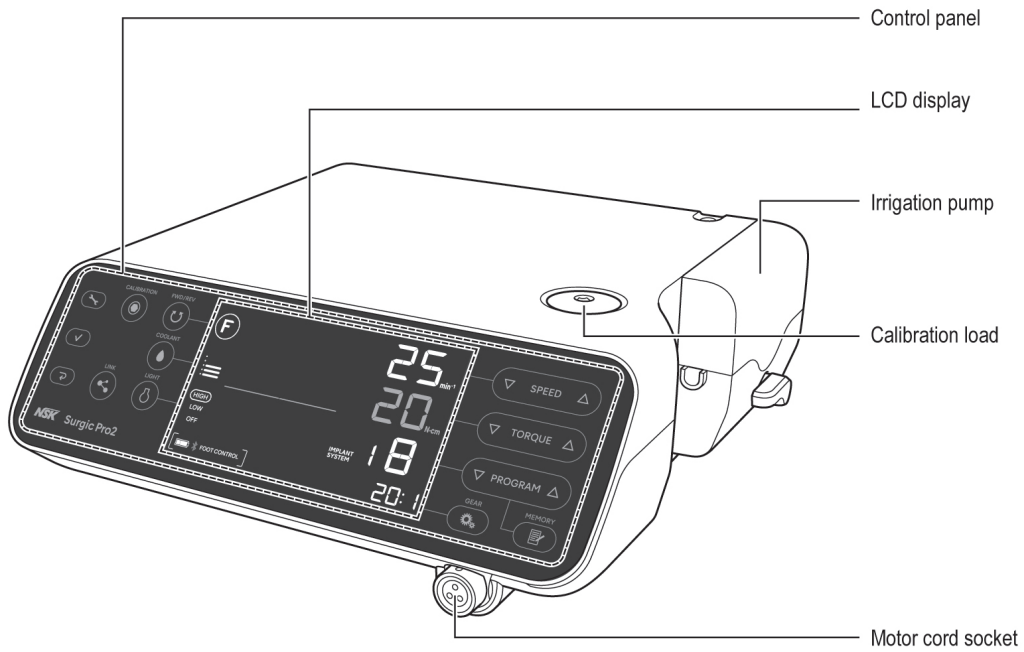
Product name
X-SG20L
SG20
X-DSG20L
X-DSG20
X-DSG20Lh
X-DSG20h

<Other products>

Product name
Link Module
Link Stand2
Osseo 100+
VarioSurg 3
VarioSurg 3 Non FT
Irrigation tube
Y-Connector


2.3 Part names

■ Control unit

2 Product Description

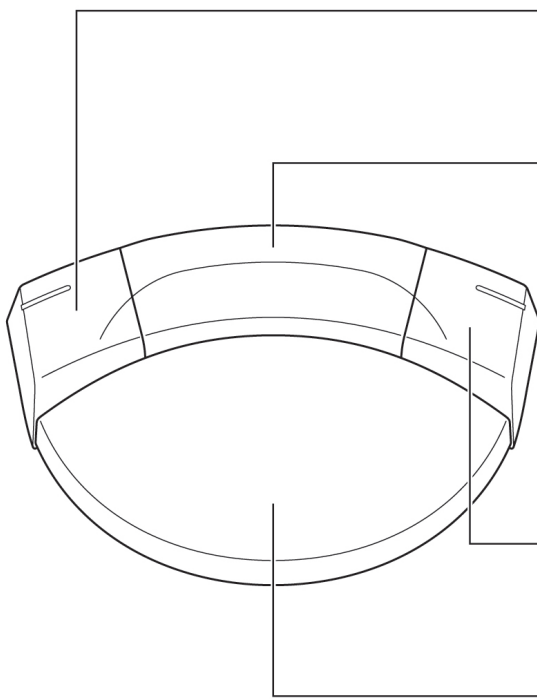■ Foot control

Left button

Select coolant solution flow level.

Press and hold for more than 2 seconds: Switching operation between this product and linked device.

Center button

Program number change

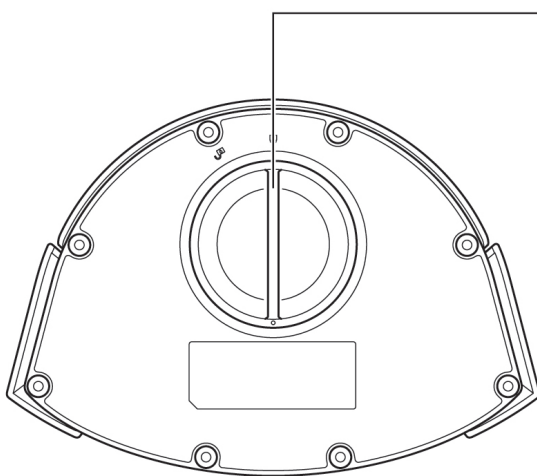
Select a program number. Depressing for about two seconds goes back one program.

Right button

Change direction of rotation

Press and hold for more than 2 seconds: Motor LED illumination while button is depressed.

Speed control pedal



Battery cover

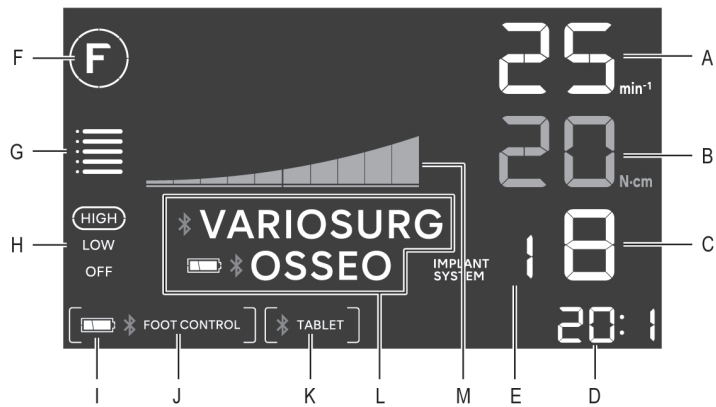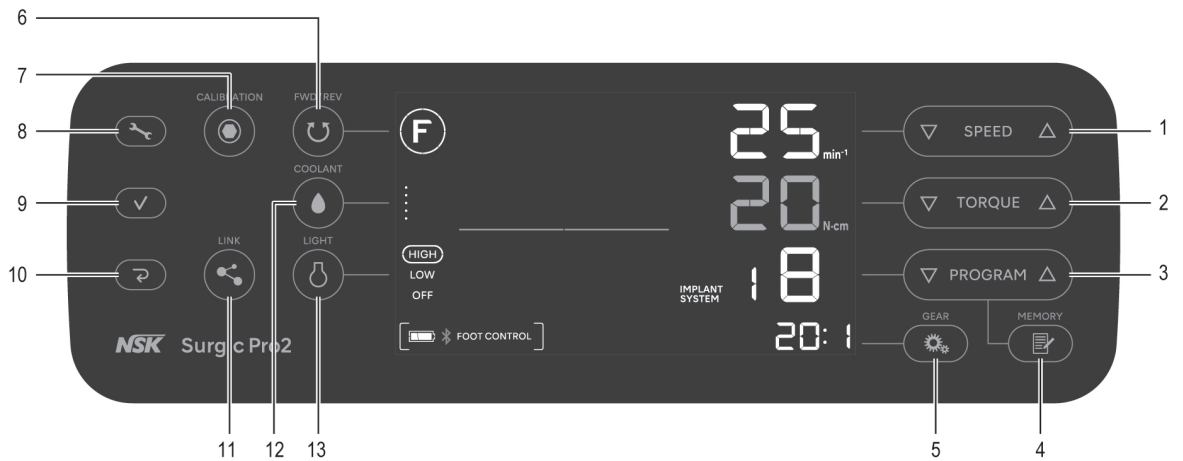NOTICE

- The functions of the buttons of the foot control noted above are factory default settings. You can change the setting as described in “7.1.2 Foot control button guidance”.
-

2.4 Control panel

No.	Key	Name	Display	Function
1		SPEED key	A	Set the motor speed.
2		TORQUE key	B	Set the torque upper limit.
3		PROGRAM key	C	Select the available programs. (1 to 8)
			E	Press ∇ and \triangle simultaneously to select the implant system. (Refer to “7.2 Implant system”)
4		MEMORY key	-	Memorize the program parameters. (gear ratio, maximum rotation speed, rotation direction, torque upper limit, coolant solution flow level, light level)
5		GEAR key	D	Select the gear ratio. (20:1, 1:1, 1:2, 1:3, 1:5)
6		FWD/REV key	F	Select the direction of rotation. \textcircled{F} : Forward rotation, \textcircled{R} : Reverse rotation
7		CALIBRATION key	-	Calibrate the handpiece. (Refer to “4.1 Calibration function”)

2 Product Description


No.	Key	Name	Display	Function
8		Setting key	-	Select the settings of the control unit, foot control or a Bluetooth® connection. (Refer to “7.1 Setting key”)
9		Enter key	-	
10		Cancel key	-	
11		LINK key	L	<p>Press and hold this key for more than 2 seconds to switch communication between VarioSurg 3 (option) and Osseo 100+(option) (only when it is connected). (Refer to “4.5 Link function”)</p> <p>* VARIOSURG : Connected VARIOSURG : Not connected, paired with VarioSurg 3 No display : Not connected</p> <p> * OSSEO : Connected OSSEO : Not connected, paired with Osseo 100+ No display : Not connected</p>
12		COOLANT key	G	Select the coolant solution flow level from 0 to 5.
13		LIGHT key	H	Select the light level. (OFF, LOW, HIGH)
-	-	-	M	<p>Displays the current operating torque. (each bar means 10%)</p> <p> : 100%  : 50%</p>
-	-	-	I	<p>Displays the battery level of the foot control. Check the battery level before and after treatment. Replace the batteries when the display of battery level is flashing.</p> <p> : 50-100%  : 10-50%  (flashing): Less than 10%  : 0%</p>
-	-	-	J	<p>Displays the status of connection with the foot control.</p> <p>[  FOOT CONTROL] : Connected</p> <p>[ FOOT CONTROL] (flashing): Not connected, paired with the foot control No display: Not connected</p>

No.	Key	Name	Display	Function
-	-	-	K	<p>Displays the status of connection between an external terminal and the control unit.</p> <p>[ TABLET] : Connected</p> <p>[TABLET] : Not connected, paired with external terminal</p> <p>No display: Not connected</p>

3 Preparation for Use


3.1 Connecting AC power cord

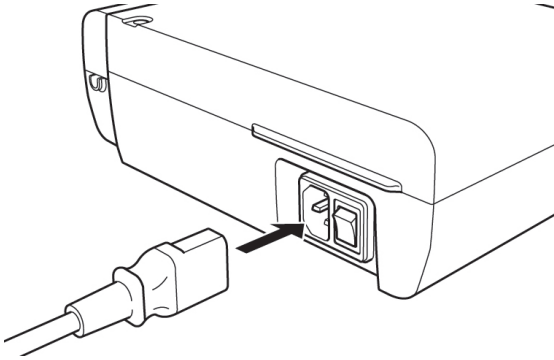Correctly align then insert the AC power cord into the AC power cord connector socket at the back of the control unit.

3.2 Installing foot control batteries

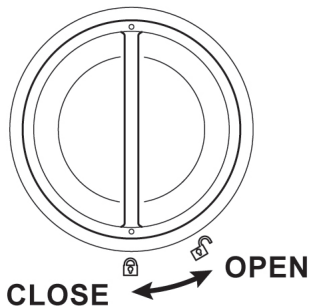1. Turn the battery cover counterclockwise.
2. Insert three AAA batteries.
3. Turn the battery cover clockwise to close.

⚠ CAUTION

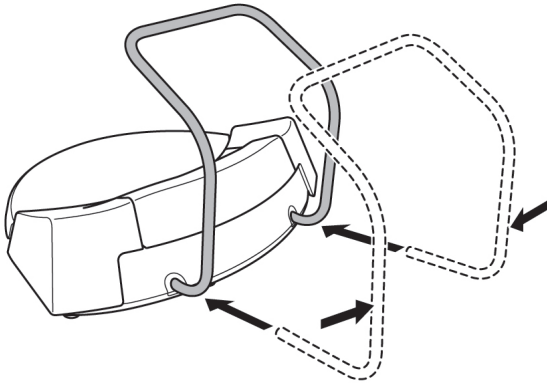
- When inserting batteries, check the marking on the battery case and make sure batteries are inserted in the correct orientation.
- Use only high-quality disposable AAA (Micro/LR03) 1.5V alkaline batteries. Using the wrong type of battery may lead to product malfunction.
- Do not mix new and old batteries or different types of batteries.
- Do not use rechargeable batteries.
- If the device will not be used for a long period of time, remove the batteries for safety reasons. Failure to do so poses a risk of product malfunction due to battery heating or leakage.
- Make sure that the O-ring attached to the cover is free of debris before closing the battery cover.

3.3 Installing foot control hanger

Squeeze the hanger and insert it aligned with the guides in the foot control holes.

Insert in the direction shown by the arrows until the hanger stops.

If removing the hanger such as when putting it in the carrying case, remove it by the reverse of the procedure for attaching. If the hanger is stiff, pull it out alternating between right and left.


⚠ CAUTION

- When holding the foot control, do not hold with the hanger retracted. It may cause the hanger to come out and the foot control to fall.


3.4 Mounting coolant solution hanger post

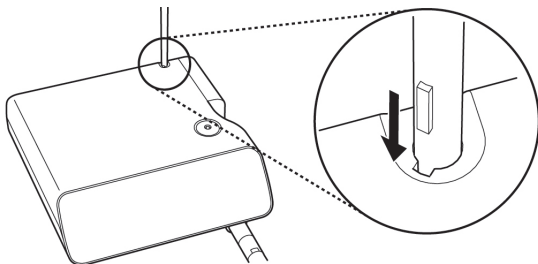Mount the coolant solution hanger post onto the holder on the control unit.


3.5 Installing irrigation tube

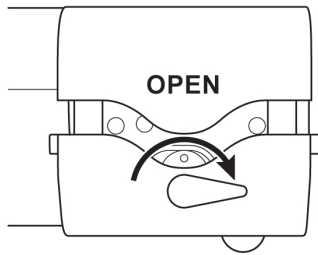
1.

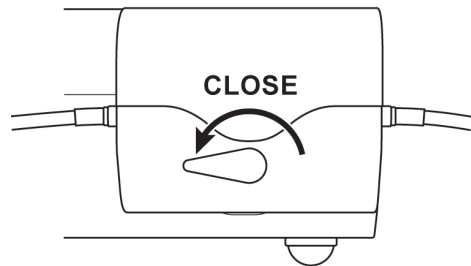3.

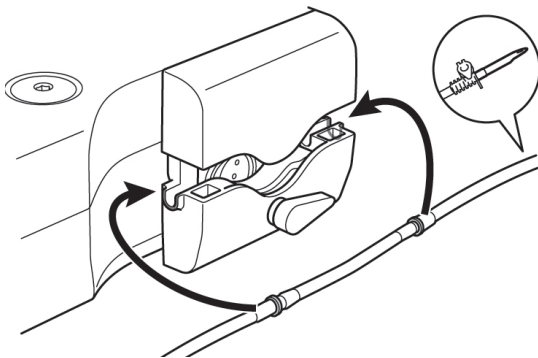
2.



⚠ CAUTION

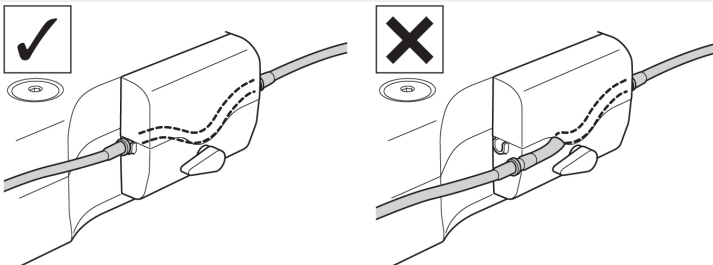- Make sure that the tube is securely set on the rollers when closing the pump cover. If the tube is not correctly positioned on the rollers and the cover is closed, the tube could be cut or sheared.
- Do not use non-NSK genuine irrigation tubing.
- When mounting the irrigation tube, be careful not to get your hand caught in the moving part of the irrigation pump.
- Do not mount the irrigation tube while the power is ON.
- Do not use the irrigation tube if the pouch is punctured or the seal is broken. It is not sterile.
- Federal law restricts this device to sale by or on the order of a physician.


3.6 Inserting irrigation tube onto and from bag / bottle

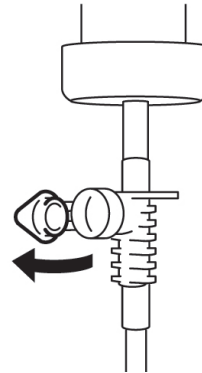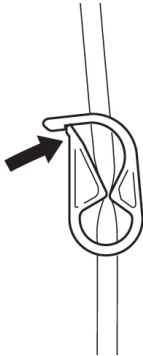1. Close the tube clamp between the irrigation tube needle and the irrigation pump.
3. Open the tube cap to supply air into the bottle.

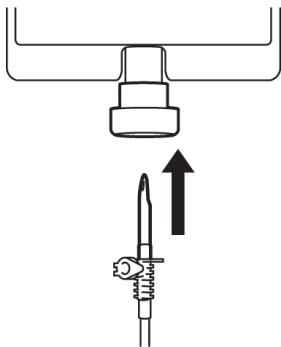2. Insert the irrigation tube needle into the bottle cap.
4. Open the tube clamp.


⚠ CAUTION

- Do not operate the irrigation pump if the tube is bent or the tube clamp is in the closed position. This could cause the tube to break or slip out of the bottle.

3.7 Connecting motor cord

Align [•] mark on the motor cord with the [•] mark on the control unit.

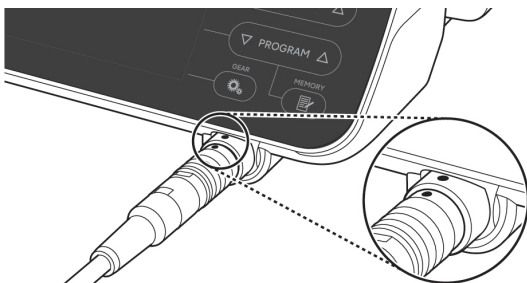Push the base of motor cord plug into the socket until the lock joint 'clicks' into place.

To disconnect the plug, pull back the lock joint and then pull further to disconnect the cord.

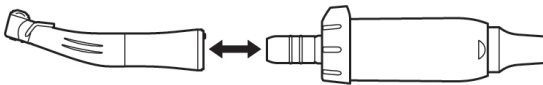
3 Preparation for Use

⚠ CAUTION

- Insert the motor cord straight and do not put excessive load on the motor cord socket as it may cause breakage or deformation.

3.8 Connecting handpiece

1. Push the handpiece on to the motor taking care to keep it in alignment.
2. Twist the handpiece until it locks into position with a clicking sound.

⚠ CAUTION

- When using the optic micromotor, use only optic handpieces. A non-optic handpiece may not precisely connect to the micromotor and it may malfunction when the micromotor is under load. Do not connect the optic handpiece to the non-optic micromotor.
- Always connect the sterilized handpiece.

NOTICE


- Refer to the handpiece manual regarding bur / drill attachment and removal and irrigation methods.

3.9 Irrigation nozzle attachment

X-SG20L and SG20 has three irrigation methods available depending on the tool and application: external, internal or both. For details on installation, refer to the included handpiece operation manual.

3.10 Attaching tube holder

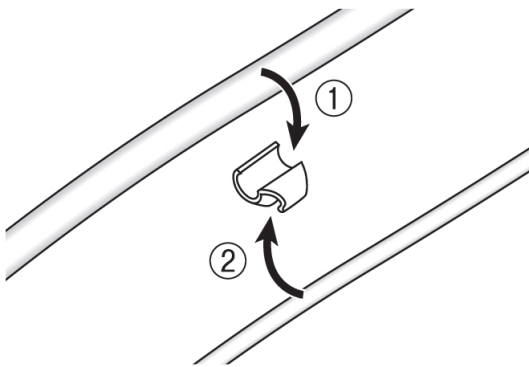Use the motor cord as a strain relief for the irrigation tube. It is easier to insert the motor cord first, then the irrigation tube.

3.11 Check before treatment

Perform the following checks before using the product on a patient to confirm that there are no abnormalities. If abnormal vibration, noise, or overheating is detected, stop using the product and contact your authorized NSK dealer.

1. Make sure each component is securely connected.
2. Make sure the control unit main power switch is in the OFF position (O side), and then connect to a commercial power socket.
3. Put the control unit main power switch in the ON position (I side).
4. Confirm compatibility as covered in "3.11.1 Compatibility check of internal spray nozzle/drill".
5. Run the motor for 1 minute and check the following.
 - Is the bur or drill free from rattling, abnormal vibration, or abnormal noise?
 - Can irrigation be done and is the coolant solution flow level appropriate?
 - Does the motor light illuminate (optic motor only)?
6. Stop the motor and make sure there is no abnormal heating of the motor and handpiece surface.

⚠ CAUTION

- Touching the control panel of an implant micromotor may increase the risk of cross-infections. In order to minimize such a risk, it is possible to use a barrier film in combination with the control unit. Do not use a film with poor operability and/or visibility.

3.11.1 Compatibility check of internal spray nozzle/drill

Internal spray nozzles included with this product do not necessarily fit all drills on the market. Follow the instructions given below for confirmation prior to use.

Failure to do so or to properly fit the internal spray nozzle into the drill may cause a saline solution leak, which will result in problems such as rust or sudden stop of equipment during use.


* To mount the drill, refer to attached handpiece operation manual.

Instructions:

1. Attach a bottle of saline solution to the control unit.
2. Connect the internal spray nozzle into the tip of the irrigation tube.
3. Insert the internal spray nozzle into the drill from the back.
4. Turn the control unit power on and irrigate at maximum flow for about 5 seconds.

Points to be checked:

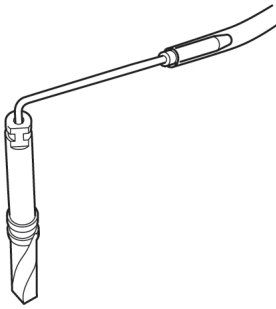- Cleanliness of the saline solution coming out from the drill: If solution is colored there could be rust inside of the drill. If so replace the drill.
- If coolant solution flow level is low, the exit point of the nozzle may become clogged with cut bone powder. Clean it or replace with a new nozzle.
- No water leakage between internal spray nozzle and drill before use: Ensure no water is leaking from the entry point of irrigation nozzle. A broken seal or no seal in the drill may be the cause. Replace the drill even if its new. Saline solution ingress into handpiece will cause malfunction.

⚠ CAUTION

- If a malfunction is detected in use, such as a saline solution leaking from the back of the contra angle head, stop using the handpiece and perform troubleshooting.

4 Operation Procedure

4.1 Calibration function


The rotation resistance of a handpiece depends on the handpiece model, condition, and internal wear of the handpiece gears. This product has a function to recognize rotational resistance when in use and make corrections so the set torque and output torque are the same.

NOTICE

- This equipment is optimized to obtain the highest calibration accuracy at a gear ratio of 20:1. For calibration USE ONLY the 20:1 reduction handpieces listed in <Contra angle handpiece> in “2.2 Applicable products”. Correct calibration cannot be done if other handpieces are used.
- Microsaw handpiece cannot be calibrated.
- Calibration should be performed only on an NSK handpiece.
- Only speed calibration is performed if the gear ratio is other than 20:1.
- Press the cancel key for 2 seconds to return to the normal screen from calibration mode.

4.1.1 Calibration preparation

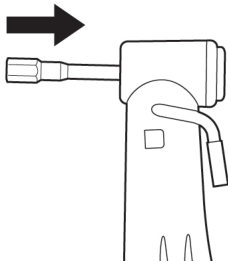1. Attach the 20:1 reduction handpiece to the motor.
2. Attach the calibration bur to the handpiece.
3. Press  (CALIBRATION) for about 2 seconds to switch to calibration mode.
A beep will sound and the display will switch to calibration mode and “CAL-TRQ” is displayed.

4.1.2 Conducting calibration

Torque calibration with no-load

1. Make sure that "L" is displayed on the LCD display.

2. Hold the motor in your hand.

3. Press  (CALIBRATION) without applying load.
4. A beep sounds and the calibration bur automatically rotates for a few seconds.
5. Calibration is complete when a beep sounds and "PASS" appears on the display. Proceed to next step.

CAUTION

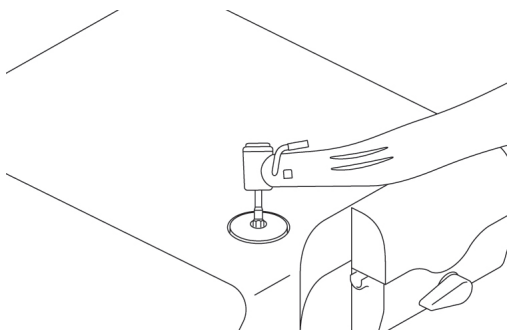
- Do not apply load during no-load calibration. If load is applied, "FAIL" will appear on the LCD display, and calibration will stop.

Torque calibration with load

1. Make sure "H" is displayed on the LCD display.

2. Plug the calibration bur into the calibration load.

3. Hold the handpiece firmly so the calibration bur is perpendicular to the calibration load.
4. Press  (CALIBRATION).
5. A beep sounds and the calibration bur automatically rotates for a few seconds.
6. Calibration is complete when a beep sounds and "PASS" appears on the display. Proceed to next step.

CAUTION

- During the process, it is important to hold the calibration bur perpendicularly without applying excessive force. Making the calibration bur tilt or pressing it may cause a loss in accuracy.

⚠ CAUTION

- During the process, be careful not to touch the control unit with the handpiece or your hands.

Speed calibration (Low speed)

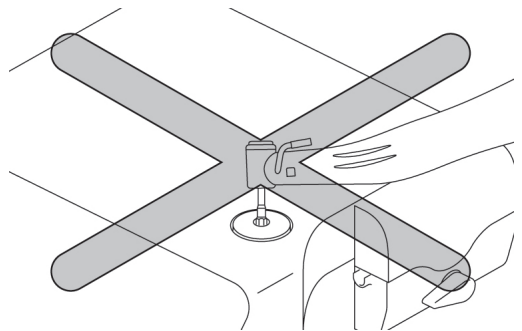
1. Remove the calibration bur from the calibration load.
2. Press  (CALIBRATION) with the bur not contacting anywhere.
3. A beep sounds and the calibration bur automatically rotates for a few seconds at low speed.
4. Calibration is complete when a beep sounds and "PASS" appears on the display. Proceed to next step.

Speed calibration (High speed)

1. Press  (CALIBRATION) with the bur not contacting anywhere.



2. A beep sounds and the calibration bur automatically rotates for a few seconds at high speed.
3. Calibration is complete when a beep sounds and "SUCCESS" appears on the display.
4. The LCD display will return to the normal screen.

Calibration fail example

- If a problem occurs such as the gear being extremely worn, correct data cannot be obtained at calibration and "FAIL" is displayed. In this case, please request repair.
- Even if a malfunction does not occur, "FAIL" is displayed if the bur comes in contact with something during no-load calibration or if you forget to apply a load at load calibration.

4.2 Starting operation

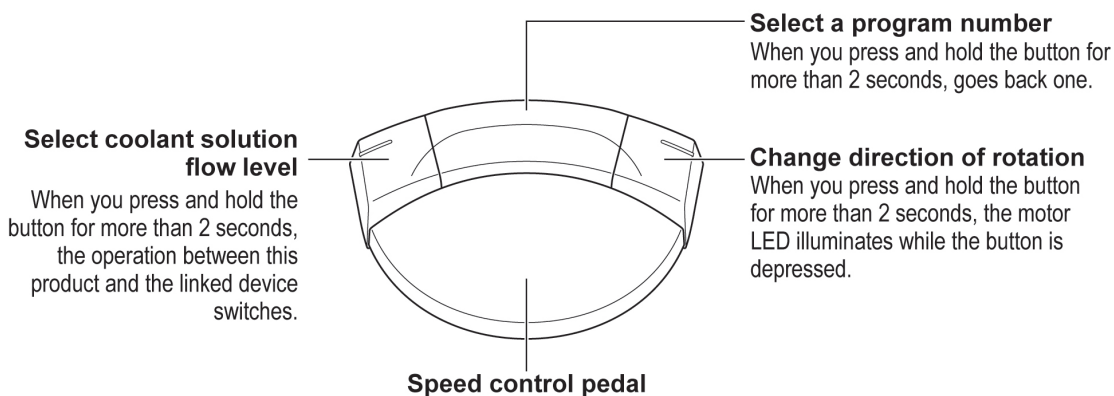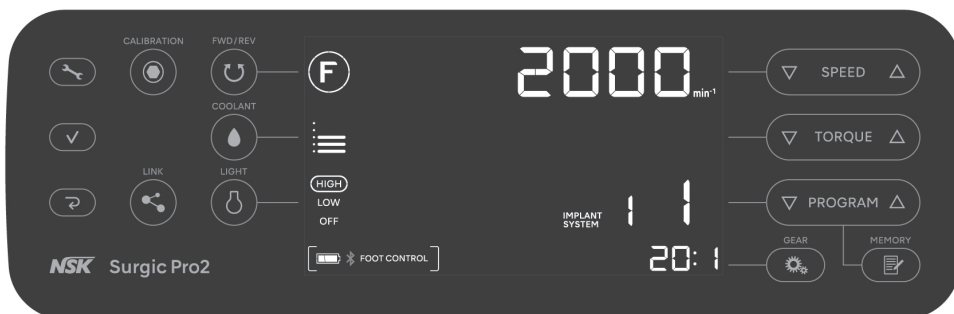Follow the instructions in “7.3 Programming system operation” and set items such as handpiece gear ratio, maximum rotation speed, direction of rotation, torque upper limit, coolant solution flow level, and light level in advance.

1. Turn on the power switch of the control unit. “CONNECT” is displayed on the LCD display as shown in the figure below, and it is automatically connected to the paired device. When the foot control is connected, the normal screen is displayed. If it is not switched to the normal screen, refer to “8.2 Problems and solutions”.
2. Press **PROGRAM** to select a program number.
3. Confirm settings such as gear ratio, maximum rotation speed, direction of rotation, torque upper limit, coolant solution flow level, and light level and change them as necessary.
4. Depress the pedal on the foot control to activate the motor.

⚠ CAUTION

- If using the motor with the light on for more than 1 minute continuously, the motor will become hot (the maximum surface temperature might reach 50 °C depending on the situation in which the system is used.) In such cases, stop using the motor until the light part cools down.
- Activation of the torque limiter: During use, when the drilling load reaches the preset torque upper limit, the integrated torque limiter automatically activates to prevent over torque. When the torque limiter activates, the motor stops after beeping and “SAFE” appears on the display. To reactivate the motor, release the speed control pedal.
- When the motor has continued operating for 15 minutes or the temperature inside the motor and the control unit rises to a certain level, “SAFE” appears on the display. To reactivate the motor, release the speed control pedal.

NOTICE

- The functions of the buttons of the foot control noted above are factory default settings. If settings are changed as described in “7.1 Setting key”, the functions will differ from those above.

4.3 Protection circuit

A protection circuit automatically functions to protect the motor and the control unit if the motor is ever overloaded. Power supply to the motor will automatically be terminated and an error code will be displayed on the control unit.

Resetting the protection circuit

To reset the protection circuit, release the speed control pedal and then depress it again.

4.4 Sleep mode

Sleep mode is a function that saves the foot control battery by shifting to sleep mode after a certain period of time when the control unit is powered on and the foot control is not operated. Refer to “7.1 Setting key” and switch FC MODE of FC SET to FC to enable sleep mode.

When the foot control enters sleep mode, “SLEEP” is displayed on the screen (LCD display). To release sleep mode, keep pressing the speed control pedal of the foot control until it returns to the normal screen. When the control unit is powered on, “PUSH FC” is displayed on the screen and press the speed control pedal of the foot control.

The time to enter sleep mode is OFF, 5 minutes, 10 minutes, 15 minutes, 20 minutes, 30 minutes, 45 minutes, and 60 minutes. Refer to “7.1 Setting key” and select the time in SLEEP of FC SET.

4.5 Link function

This product can be connected using Bluetooth to NSK VarioSurg 3 connected to the Link Module (option), Osseo 100+(option), or an external terminal such as iPad.

1. Press and hold for 2 seconds or more to switch to setting mode, and “BLE SET” is displayed.
2. Press .
3. Press to select sub menu.

Foot control

Refer to “7.1.1 Pairing with foot control”.

Osseo 100+

Refer to “4.5.1 Osseo integration monitoring device (Osseo 100+)”.

VarioSurg 3

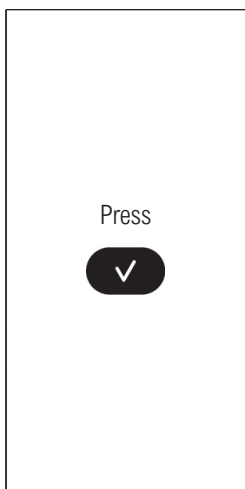
Refer to “4.5.2 Link Module”.

External terminal

Refer to “4.5.3 External terminal”.

NOTICE

- Refer to the operation manual of the device to be connected for the use method and use environment.
- If pairing is not completed within 5 minutes, “NG” is displayed.
If this happens, repeat pairing operation again from step 1.

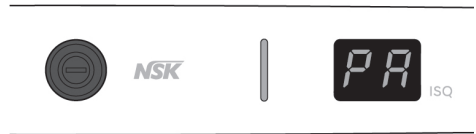
4.5.1 Osseo integration monitoring device (Osseo 100+)

The Osseo integration Monitoring Device, “Osseo 100+” measures the implant stability quotient (ISQ) contactlessly and can be connected using Bluetooth, permitting ISQ confirmation over a shared interface.

The ISQ score measured by Osseo 100+ can be displayed on the control unit of this product.

Refer to the Osseo 100+ manual for the method of use.

1. Press the operating key to turn on Osseo 100+.
3. Confirm the “PA” is displayed.

2. Press and hold the operating key at least 3 seconds.

4. Pairing is successful when this product beeps and displays “OK”.

Make sure  icon below is displayed and the device is connected.

5. Press  twice to end setting mode.

⚠ CAUTION

- If any abnormalities are found with Osseo 100+, refer to the Osseo 100+ operation manual.
- Do not use the ISQ value displayed on the screen of this product for diagnostic purposes. The value is displayed for information only.

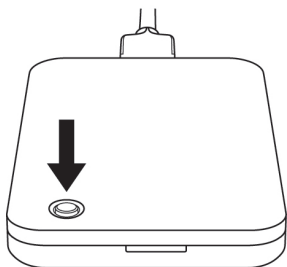
4.5.2 Link Module

The link function enables two systems (NSK VarioSurg 3 and this product) to be controlled by the foot control. Connecting Link Module to the VarioSurg 3 enables the function to link with this device to be used.


Refer to the OM-E0859E Link Module OPERATION MANUAL listed in “1.2 Related manuals” for the method of use.

1. Turn on the VarioSurg 3.
2. Press and hold the button on Link Module.
4. Pairing is successful when this product beeps and displays “OK”.

Make sure  icon below is displayed and the device is connected. Also the LED of the Link Module turns green.

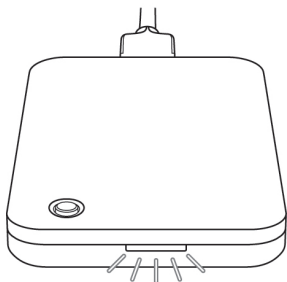3. When the LED turns blue, release the button.

5. Press  twice to end setting mode.

CAUTION

- If any abnormalities are found with VarioSurg 3 and Link Module, refer to the respective operation manuals.
- If the button is not released after the LED turns blue, the LED blinks white and pairing stops. To pair, release the button, and after the LED lights white, proceed from step 2.

NOTICE

- Connect the Link Module to the VarioSurg 3 beforehand.

4.5.3 External terminal

By connecting the control unit to an external terminal (such as iPad) with a dedicated application installed, detailed procedural data such as rotation speed and torque can be displayed in real time.

Procedural data can also be stored on the external terminal in PDF or CSV files.

Refer to the operation manual in the application for how to use the application.

Download the dedicated application from the following.

www.nsk-dental.com/qr/app-surgicpro2

Refer to OM-DE1112MA for compatible terminals.

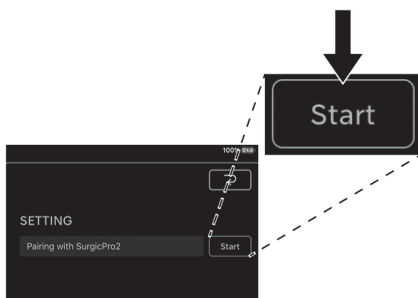
1. Launch the App.


2. Touch the icon.

3. Touch "Start".

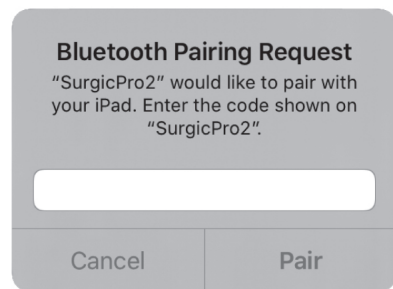4. Enter the code "123456" and touch "Pair".

5. Pairing is successful when this product beeps and displays "OK".
Make sure  icon below is displayed and the device is connected.

6. Press  twice to end setting mode.

⚠ CAUTION

- The data saved on the external terminal is for reference of operation history, do not use it for diagnostic purposes.

5 Post-use Maintenance

5.1 Preparation

1. Wear eye protection, a mask, and gloves to prevent infection.
2. Turn off the power switch of the control unit.
3. Remove the bur/drill from the handpiece.
4. Remove the handpiece from the motor.
5. Remove the motor cord from the control unit.

⚠ WARNING

- For cleaning and disinfection, use disinfectant ethanol or disinfectant isopropanol that does not contain any additives (hereinafter referred to as disinfectant alcohol). Use of disinfectants other than disinfectant alcohol may cause discoloration or cracking.
- Failure to properly maintain the products may cause infection, product failure, or overheating and injuries. Be sure to clean and disinfect the products immediately (within 1 hour) after use to remove residue.

⚠ CAUTION

- Use appropriate disinfectant alcohol according to the regulations and guidelines of each country.
- Do not clean, immerse, or wipe with electrolyzed-oxidizing water (strong acidic water, super acidic water), strong acidic and strong alkaline liquid solution, solvent containing chlorine, benzene or thinner.
- Follow any additional local directives, standards, and guidelines for cleaning, disinfection and sterilization.
- For maintenance of the handpiece, internal spray nozzle, nozzle holder and cleaning wire, refer to the OM-C0595E/ OM-DC0929MA manual listed in “1.2 Related manuals”.
- The product lifetime and product's ability to operate correctly are mainly determined by mechanical stresses during use and the chemical influences due to processing. If you detect any abnormality such as center run-out of bur/drill, heat generation, vibration or noise, or if the exterior markings (serial number or model number) have become unreadable before or during use, stop using the product immediately and request your nearest NSK dealer, or one of the NSK branches for a periodic inspection*.

Disappearance of exterior markings is not covered under warranty for repair.

When the product has been reprocessed 250 times or has been in use for one year, it is recommended that a periodic inspection be performed by an authorized NSK service*.

*Refer to the link below:

5.2 Cleaning and disinfection

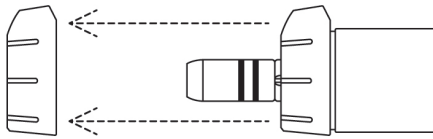
Manual Method

- For the motor with motor cord and handpiece stand


Cleaning

1. Remove the tube holder (gray) from the motor with motor cord.

For cleaning and disinfection of the tube holder (gray), refer to “■ For the tube holders, tube holder (gray) and protection plug”.

2. Attach the protection plug to the motor with motor cord.

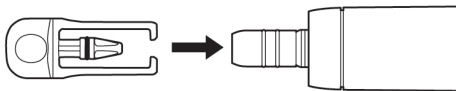3. Clean the surface of the products under running water with a soft bristled brush under the following conditions.

Time: 1 minute for the motor with motor cord, 15 seconds for the handpiece stand

Water temperature: ≤ 38 °C

Water quality: Drinking water

Water flow rate: ≥ 3.5 L/min

4. Remove any remaining moisture from the product using a dry cloth or filtered clean and dry compressed air (≤ 0.35 MPa) ensuring it is dry internally and externally.
5. Remove the protection plug from the motor with motor cord. Under appropriate lighting (≥ 500 lx), inspect the products for contamination. If any visible contamination remains, attach the protection plug to the motor with the motor cord again and repeat the process until it is visually clean following step 3.

Disinfection

1. Wipe the surface of the products with disinfectant wipes (MinutenWipes from ALPRO) for 2 minutes each.

- For the tube holders, tube holder (gray) and protection plug

Cleaning

1. Clean the surface of the products under running water with a soft bristled brush for 15 seconds.

Water temperature: ≤ 38 °C

Water quality: Drinking water

Water flow rate: ≥ 3.5 L/min

2. Proceed to the section “5.3 Sterilization”.

2. Remove any remaining moisture from the product using a dry cloth or filtered clean and dry compressed air (≤ 0.35 MPa) ensuring it is dry internally and externally.
3. Under appropriate lighting (≥ 500 lx), inspect the products for contamination. If any visible contamination remains, repeat the process until it is visually clean following step 1.

5 Post-use Maintenance

Disinfection

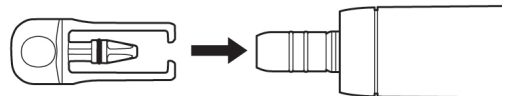
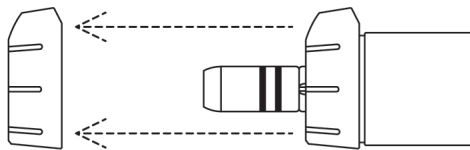
1. Spray the WL-cid (ALPRO) on to the products. Wait for the WL-cid (ALPRO) to take effect following the instructions provided by the manufacturer.
 2. Proceed to the section “5.3 Sterilization”.
- For control unit and foot control
1. Wipe the surface of the products with a cloth moistened with disinfectant alcohol.
 2. Store the products in a place where they can be kept clean and dry until next use.

Automated cleaning and disinfection

- For motor with motor cord, tube holders, tube holder (gray), calibration bur, protection plug and handpiece stand

 NSK products with this symbol are compatible with washer-disinfector.

1. Remove the tube holder (gray) from the motor with motor cord.
2. Attach the protection plug to the motor with motor cord.

CAUTION

- Always use the Protection Plug when using washer-disinfector.

3. Place the products into the washer-disinfector.
4. Run the washer-disinfector cycle to clean and disinfect the products. Use a Miele washer-disinfector PG8581 (or equivalent device).
5. Remove the products from the washer-disinfector.
6. Remove any remaining moisture from the product using a dry cloth or filtered clean and dry compressed air (≤ 0.35 MPa) ensuring it is dry internally and externally.
7. Under appropriate lighting (≥ 500 lx), inspect the products for contamination. If any visible contamination remains, repeat the process until it is visually clean following step 3.
8. Remove the protection plug from the motor with motor cord. Proceed to the section “5.3 Sterilization”.

Due to the variation in cleaning/disinfection cycles and cleaning agents available from the various manufacturers, NSK has validated the Miele washer-disinfector PG8581 Cycle Vario TD (Refer to the table below) for automatic processing/reprocessing of the products.

Only validated procedures must be used for cleaning and sterilization. When using reprocessing procedures that are different from those described in this manual, those procedures must be validated by the respective practice or hospital using washer-disinfector that conforms to EN ISO 15883-1.

Step1	Precleaning with cold tap water for 1 minute.
Step2	Cleaning for 5 minutes using demineralized water mixed with 0.5% cleaner at 55 °C.

Step3	Rinsing with demineralized water for 1 minute.
Step4	Thermodisinfection with demineralized water at 93 °C for 5 minutes.
Step5	Dry the products for 5 minutes at 80-100 °C.

⚠ CAUTION

- Products must be removed from the washer-disinfector straight away (within 1 hour) after the cleaning, disinfecting and drying cycle is complete to prevent corrosion.
- After using the washer-disinfector, dry the product completely, because remaining moisture may cause internal corrosion, etc.
- After cleaning and disinfection of the products, be sure to perform “5.3 Sterilization”.
- Use a cleaning and disinfection device (washer-disinfector) complying with EN ISO 15883-1 (e.g. Miele washer-disinfector PG8581), that operates with a maximum pH value of 10.5 (e.g. neodisher MediClean, Dr. Weigert). Proof of suitability for the process must be obtained from the washer-disinfector manufacturer.
- Automated cleaning and disinfection should be done in accordance with the operation manual of the washer-disinfector.

5.3 Sterilization

After each patient treatment, sterilize the products as follows.

Sterilizable equipment:

Handpiece, motor with motor cord, handpiece stand, tube holders, tube holder (gray), protection plug, and calibration bur.

1. Insert the products into a sterilization pouch that conforms to EN ISO 11607-1 and seal the pouch.
2. Sterilize using the following conditions.

	Gravity displacement	Pre-vacuum (dynamic air removal)
Temperature	132 °C	134 °C
Full cycle time	15 min	3-18 min

⚠ CAUTION

- Follow local rules, regulations, and guidelines regarding the reprocessing of devices.
- Do not use sterilization pouches containing water-soluble adhesive components such as PVA (polyvinyl alcohol). Adhesive components washed out during sterilization may enter the product, causing problems such as poor rotation, sticking and failure to operate correctly.
- Do not autoclave the product with other instruments even when it is in a pouch. This is to prevent possible discoloration and damage to the product from chemical residue on other instruments.
- Be sure to use sterilizers that can perform sterilization up to 135 °C. In some sterilizers, the chamber temperature may exceed 135 °C. Do not use these sterilizers as failure of the product could occur. Contact the sterilizer manufacturer for detailed information about cycle temperatures.
- Do not heat or cool the product too quickly. Rapid change in temperature will cause damage to the product.
- The irrigation tube is for single use product and cannot be autoclaved.
- Do not touch the product immediately after autoclaving as it will be very hot.
- Steam sterilization is recommended for the product. The validity of other sterilization methods (such as plasma sterilization or EOG sterilization) is not confirmed.

5 Post-use Maintenance

CAUTION

- Immediately after sterilization is complete (within 1 hour), remove the products from the sterilizer. Failure to do so may cause corrosion.
-

NOTICE

- NSK recommends use of a Class B sterilizer as stated in EN 13060.
 - For details, confirm the manufacturer's sterilizer operation manuals.
-

5.4 Storage

Store the product in a place where it is kept clean and keep it in a sterilization pouch until it is used next.

CAUTION

- Keep the product in suitable atmospheric pressure, temperature, humidity, ventilation, and sunlight. The air should be free from dust, salt and sulphur.
 - Store the product in a well ventilated place out of direct sunlight and within the range of temperature, humidity and pressure specified in "10 Specification".
 - Sterilization is not guaranteed after the sterilization retention period specified by the manufacturer and seller of the sterilization pouch has elapsed. If the sterilization retention period has elapsed, perform sterilization again with a new sterilization pouch.
-

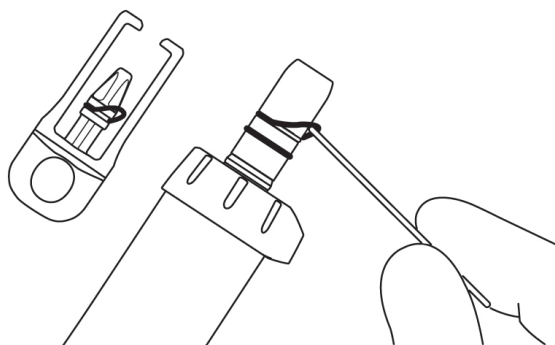
6 Maintenance

6.1 Replace O-ring

If the O-ring has become worn or damaged, the handpiece may become difficult to connect or if the handpiece is rattling, replace the O-ring.

Remove the O-ring with a pointed tool, and then mount the new O-ring into the groove.

Also replace the O-ring of the protection plug as well.

6.2 Periodical maintenance checks

Every 3 months perform periodical maintenance checks, referring to the check sheet below. If any abnormalities are found, contact an authorized NSK dealer.

Points to check	Details
Control unit operation	There should be no abnormal mechanical sound. LCD display is displayed appropriately. The control panel does not operate abnormally.
Foot control operation	Check the battery level displayed on the control unit LCD display. If battery level is low, replace the batteries. Fully depress the speed control pedal and confirm that the motor is rotating at the set speed. Depress each button and make sure that function corresponding to each button can be confirmed. If the foot control does not function properly, perform calibration. (Refer to “6.3 Foot control calibration”) Remove the foot control battery cover and make sure its O-ring is not damaged or worn. If damage or wear is seen, contact an authorized NSK dealer for replacement.
Foot control hanger	Make sure the foot control hanger does not easily come off the foot control.
Coolant solution flow	Coolant solution flows appropriately with no leakage.
Motor with motor cord	Operate for 1 minute and make sure there is no abnormal heating, sound, or vibration. Make sure the motor LED illuminates normally.
Handpiece	Inspect the handpiece according to its manual.

Additional safety information

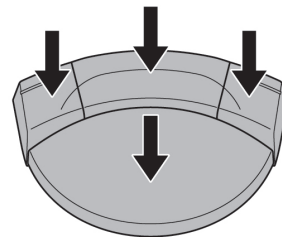
Periodical safety testing shall be performed every two years following IEC 62353 (EN 62353), in accordance with test procedure OM-DE0938EN. As periodical safety testing is limited exclusively to professional dental installers who have been trained and qualified under the control of NSK, contact your authorized NSK dealer.

6.3 Foot control calibration

Foot control pedal and button response may become poor due to wear of parts over time. In that case, perform calibration as described below.


1. Press and hold for 2 seconds or more to switch to setting mode.
2. Press to select "FC SET".
3. Press .
4. Press to select "FC CAL".
5. Press to select "ON".
6. Press and confirm that "---:--" is displayed.
7. Fully depress the pedal and all of the buttons on the foot control.

8. When "000:0" is displayed, calibration is successful. Press .
9. Make sure that "OK" is displayed on the LCD display.
10. Press twice to end setting mode.

⚠ CAUTION

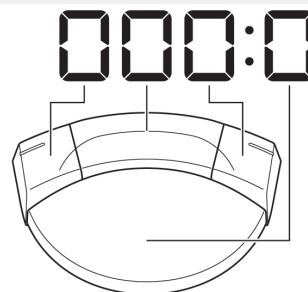- Do not press the enter key when other than "000:0" is displayed.
If calibration is not performed correctly, "NG" is displayed on the LCD display.
In this case, start again from step 1.

7 Setting

7.1 Setting key

Control unit and foot control operational settings can be changed.

For setting items and factory settings, refer to “<Setting value>”.

- Press and hold for 2 seconds or more to switch to setting mode, and “BLE SET” is displayed.
- Press to select main menu.
BLE SET : Settings for Bluetooth connection with control unit
UNIT SET : Settings for control unit
FC SET : Settings for foot control
- Press .
- Press to select sub menu. Refer to “<Setting value>”.
- Press to change the setting value.
- Press to finalize the setting value. To change other sub menu setting items, repeat procedure from step 4 to 6.
- Press to return to main menu selection screen. To change other main menu setting items, repeat procedure from step 2 to 6.
- Press to end setting mode.

<Setting value>

Main menu	Sub menu	Setting details	Setting range	Factory setting
BLE SET	PAIRING	Pair the foot control with the control unit	FC	-
		Pair Osseo 100+ with the control unit	OS	-
		Pair the Link Module with the control unit	LK	-
		Pair an external terminal with the control unit	PC	-

Main menu	Sub menu	Setting details	Setting range	Factory setting		
BLE SET	BLE OFF	Unpair the foot control from the control unit	FC	-		
		Unpair Osseo 100+ from the control unit	OS	-		
		Unpair the Link Module from the control unit	LK	-		
		Unpair the external terminal from the control unit	PC	-		
UNT SET	BZ VOL	Volume for beep that sounds at key operations and when activating.	H: high L: low	H		
	LCD BL	Operation panel brightness	1 - 10 Larger value: Brighter	10		
	TS LV	Touch sensor sensitivity level (When sensitivity is poor such as when using barrier film or surgical gloves)	1 - 3 Larger value: More sensitive	1		
	CL LV1	Coolant solution flow level (per display)		1 - 100%	20	
	CL LV2				1 - 100%	40
	CL LV3				1 - 100%	60
	CL LV4				1 - 100%	80
	CL LV5				1 - 100%	100
F-RESET	Initializing settings When display switches from off to on, press the enter key to restore settings to factory settings.					
FC SET	BUT ST	Foot control button guidance (Refer to “7.1.2 Foot control button guidance”)	Refer to notes on the left			

7 Setting

Main menu	Sub menu	Setting details	Setting range	Factory setting
FC SET	PDL-SET	Set the motor rotation specifications according to the amount the foot control speed control pedal is depressed. VAR: Rotation speed changes within the set rotation speed range, according to the amount the pedal is depressed FIX: When the pedal is depressed more than 50 %, the motor rotates at the set speed	VAR, FIX	VAR
	FC CAL	Calibrate the foot control pedal and buttons. (Refer to “6.3 Foot control calibration”)		-
	FC MODE	AUT: Set the foot control to the automatic connection. FC: Enable sleep mode. (It may be necessary to keep pressing the speed control pedal of the foot control to connect the control unit and the foot control automatically.) (Refer to “4.4 Sleep mode”)	AUT, FC	AUT
	SLEEP	Set the time until the foot control shifts to sleep mode. *This can be set only when FC MODE is changed from AUT to FC. (Refer to “4.4 Sleep mode”)	OFF 5m 10m 15m 20m 30m 45m 60m	5m

CAUTION

- When setting values are changed, press the enter key to finalize. If you press the SPEED key, TORQUE key, or Cancel key without finalizing with the Enter key or if you turn the power off, setting details will not be applied.

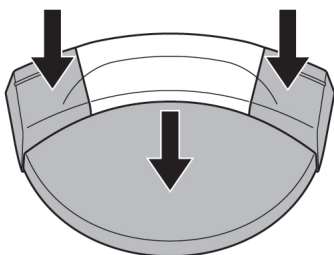
NOTICE

- Initializing settings returns “UNT SET” in setting mode and program settings to factory default. “BLE SET” and “FC SET” in setting mode are not initialized. Make a memo or other record of the current settings as necessary.
- Finalized setting details will be retained even if the power is turned off.

7.1.1 Pairing with foot control

1. Perform the operation procedure of “4.5 Link function”.
2. Press .
3. Press and hold the two buttons and the pedal of the foot control shown below simultaneously for approximately 10 seconds.
4. Pairing is successful when a beep sounds and “OK” is displayed on the LCD display. Make sure battery and  icons below are displayed, and the device is connected.

5. Press  twice to end setting mode.

⚠ CAUTION

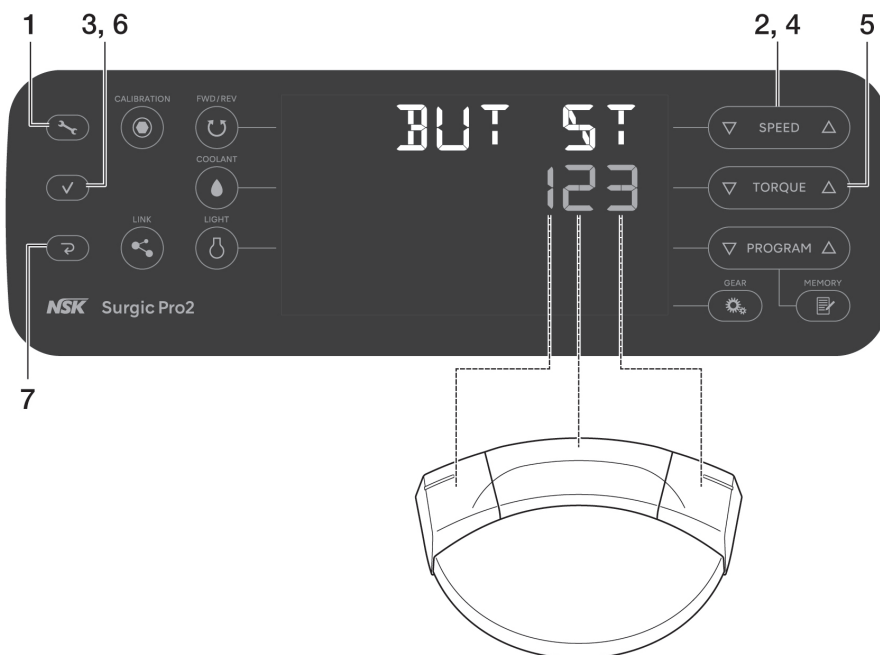
- If pairing cannot be accomplished, replace the batteries with new ones.
- Pairing may not be possible if something is obstructing the signal between the foot control and control unit.

NOTICE

- Pairing can be done within 10 minutes after inserting the battery.

7 Setting

7.1.2 Foot control button guidance

1. Press and hold  for 2 seconds or more to switch to setting mode.
2. Press  to select "FC SET".
3. Press .
4. Press  to select "BUT ST".
5. Press  to change setting value. Refer to "<Setting values of foot control button>".
6. Press .
Press  each time to select the buttons to which functions are assigned in the order of Left button → Center button → Right button → Left button ...
7. Press  twice to end setting mode.

<Factory settings>

Left button	Center button	Right button
1	2	3

<Setting values of foot control button>

Setting value	Function	
	Press and release button	Press and hold button more than 2 seconds
0	Invalid	Invalid
1	Coolant solution flow level selection (up)	Switching operation between this product and linked device
2	Program number selection (up)	Program number selection (down)
3	Direction of rotation switching	Motor LED illumination while button is depressed (Motor will not rotate)
4	Rotation speed setting (up)	Rotation speed setting (down)
5	Rotation speed setting (up)	Rotation speed setting (fast forward to up-direction)
6	Rotation speed setting (down)	Rotation speed setting (fast forward to down-direction)
7	Torque setting (up)	Torque setting (down)

Setting value	Function	
	Press and release button	Press and hold button more than 2 seconds
8	Torque setting (up)	Torque setting (fast forward to up-direction)
9	Torque setting (down)	Torque setting (fast forward to down-direction)
A	Coolant solution flow level selection (up)	Coolant solution flow level selection (0: no coolant solution flow)
B	Program number selection (up)	Implant system selection (up)
C	Implant system selection (up)	Implant system selection (down)
D	Motor LED brightness selection	
E	Motor LED illumination while button is depressed (Motor will not rotate)	
F	Irrigation while button is depressed (Motor will not rotate)	

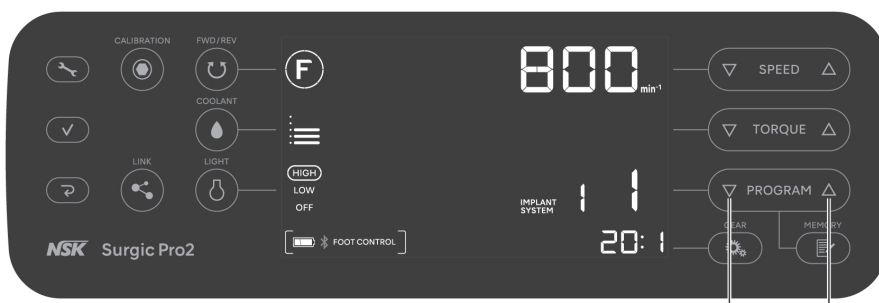
7.2 Implant system

Eight programs can be set for each of the eight different implant systems.

When you set the program, refer to values the implant manufacturer recommends.

How to change implant system

Press the ∇ and Δ on ∇ PROGRAM Δ simultaneously to change the implant system.

* Factory settings: Refer to “10.5 Factory settings (implant system)”.

7.3 Programming system operation

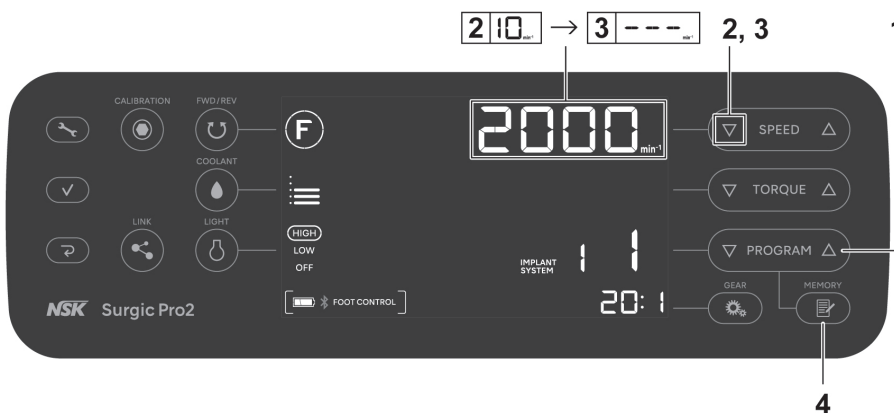
Eight programs can be memorized and recalled according to treatment procedures.

1. Press **PROGRAM** to select a program number.
2. Press **FWD/REV** to select the rotation direction:
 - F** Forward rotation (clockwise)
 - R** Reverse rotation (counterclockwise)
 - * Beep sounds constantly when reverse rotation is selected.
3. Press **GEAR** to select gear ratio for handpiece used.
4. Press **SPEED** to set the maximum rotation speed.
 - * When rotating: displays actual speed.
 - * When stopped: displays maximum speed set.
 - * Beep sounds when upper or lower limit is reached.
5. Press **TORQUE** to set torque upper limit (20:1 Reduction contra with rotation speed of 100min⁻¹ or less only).
 - * Stops once at 50N•cm. If setting to greater than 55N•cm, press **TORQUE** again.
 - * Beep sounds when upper or lower limit is reached.
 - * No display when the gear ratio is other than 20:1.
6. Press **COOLANT** to select coolant solution flow level.
7. Press **LIGHT** to select light brightness.
8. Press **MEMORY** to memorize setting (Memorizing is complete when long beep sounds).

7.3.1 Program skip function

When switching program numbers with the foot control, numbers not to be used can be memorized in advance and those selected to be skipped.

1. Press **PROGRAM** to choose program numbers that you want to skip.
2. Press **SPEED** to set the minimum rotation speed.
3. Press **SPEED** again to display “---” on the LCD display.
4. Press and hold **MEMORY** more than 1 second until it beeps to memorize.

Canceling the program skip function

1. Press **PROGRAM** to choose program numbers that you want to cancel the skip function.
2. Press **SPEED** to set a speed other than “---”.
3. Press and hold **MEMORY** for more than 1 second until it beeps to memorize.

8 Troubleshooting

8.1 Error code

If an error occurs and the motor stops, an error code is displayed on the LCD and a beep sounds to make it easy to check the status of the control unit and identify the cause of the error. The error code will blink until the error is alleviated.

Depress the speed control pedal of foot control or press the Cancel key to perform an error check. If no error is found, the error condition can be canceled. If the foot control is not connected, press the Cancel key to cancel the error condition.

Error

code	Type of error	Cause of error	Action
E0	System error	Self-check error.	Request repair.
E1	Motor voltage/ current detection error	Extended use under heavy load. Short circuit in the motor. Motor cord (power line) failure.	Electrical contact may be insufficient. Securely re-connect the motor cord. If the error cannot be eliminated, request repair.
E2	Power supply error	Insufficient contact of the power cord. Failure of the power supply circuit.	Securely re-connect the power cord. If the error cannot be eliminated, request repair.
E3	Motor sensor error	Motor sensor failure. Unplugged motor cord. Motor cord (signal line) failure. Ingress of water into the motor.	Electrical contact may be insufficient. Securely re-connect the motor cord. If the error cannot be eliminated, request repair. (Make sure to attach a protection plug when cleaning in a washer-disinfector.)
E4	Control unit interior overheating error	Overheating by extended use under heavy load. Use immediately after the control unit is placed at high temperature (such as inside a car under the blazing sun or in a locker exposed to direct sunlight).	Allow it to cool down before use. To allow heat to be sufficiently radiated, the periphery of the control unit should be well-ventilated wherever possible. If the error cannot be eliminated, request repair.
E5	Breaking error	Abnormal voltage generated in the start / stop switch circuit. Failure in the start / stop switch circuit.	When rotation and stop are repeated in a short period of time, a circuit to suppress the sudden acceleration at the start of rotation may be activated. Wait for a few seconds before use. If the error cannot be eliminated, request repair.
E6	Motor rotation failure error	Bur or drill attachment failure. Handpiece failure. Motor failure.	The chuck may be opened, or may not be sufficiently closed. Securely close the chuck. If the error cannot be eliminated, request repair.
E7	Irrigation pump error	Incorrectly fitted irrigation tube caught on pump roller. Irrigation pump failure.	Check the irrigation tube. If the error cannot be eliminated, request repair.

Error code	Type of error	Cause of error	Action
E8	Bluetooth communication error	Failure of the Bluetooth communication circuit.	Request repair.
E9	Foot control error	Failure of the foot control sensor (Hall IC). Self-check error.	Request repair.
E10	Motor LED error (Optic motor)	Adhesion of the saline solution to the motor LED. Failure of LED circuit.	Remove the saline solution completely. If the error cannot be eliminated, request repair.
E11	Control panel error	Failure of the LCD display driver. Failure of the electrostatic sheet.	Request repair.
E12	Link Module error	Failure of the link cable. Failure of the Link Module circuit. Failure of Bluetooth communication circuit in the Link Module.	Turn off the VarioSurg 3 and re-connect the Link Module securely. If the error cannot be eliminated, request repair.
E13	Osseo 100+ error	It is affected by electromagnetic interference from electronics Influence of electromagnetic interference waves emitted by electronic equipment. Failure of Bluetooth communication circuit in Osseo 100+.	Keep the sufficient distance between the tip of Osseo 100+ and the electronic equipment. If the error cannot be eliminated, request repair.
E14	External terminal error	OS of external terminal is not supported. Failure of the external terminal.	Check the supported OS in the application. Also, reinstall the application. If the error cannot be eliminated, request repair.

8.2 Problems and solutions

When a problem is detected, check the following again before requesting a repair. If none of these is applicable or if the trouble is not remedied even after an action has been taken, a failure of this product is suspected. Contact an authorized NSK dealer.

Problem	Cause of problem	Action
Screen is blank.	Plug improperly inserted to power source	Check to make sure the plug is in.
	Malfunction of power switch	Request repair.
The maximum rotation speed and "--" are alternately displayed on the speed display.	Power turned on while pressing the foot control	Remove foot from foot control.
	Foot control malfunction	If the problem is solved by removing the batteries, it is a foot control malfunction. If not, it is a control unit malfunction. Please contact an authorized NSK dealer.
	Control unit malfunction	

8 Troubleshooting

Problem	Cause of problem	Action
After the power is ON, "CONNECT" or "PUSH FC" is displayed and the screen is not switched.	Foot control is not connected.	<p>Keep pressing the speed control pedal of the foot control while "PUSH FC" is displayed on the LCD display. If the screen does not change after 5 minutes, press and hold the Cancel key. The normal screen is displayed.</p> <p>Refer to "3.2 Installing foot control batteries" and replace the foot control batteries to check if it can be connected.</p> <p>If it cannot be connected after replacing the batteries, Refer to "7.1.1 Pairing with foot control" for pairing. If this does not resolve the issue, please contact an authorized NSK dealer.</p>
Will not rotate even if the foot control is pressed.	Foot control is not connected.	<p>Perform the following after confirming that nothing is blocking the signal between the foot control and control unit.</p> <p>Turn the main unit power off and on again to connect to the control unit.</p> <p>Refer to "3.2 Installing foot control batteries" and replace the foot control batteries with new ones.</p> <p>Refer to "7.1.1 Pairing with foot control" and perform pairing.</p> <p>If this does not resolve the issue, please contact an authorized NSK dealer.</p>
	Foot control pedal and button response is poor due to not being used for a long period of time.	Refer to "6.3 Foot control calibration" and calibrate the foot control.
	Malfunction of control unit or foot control	If the problem is not solved by checking foot control batteries or foot control connection status, system repair is required.

Problem	Cause of problem	Action
"COM ERR FC" is displayed.	Foot control not connected	<p>Press the Cancel key to cancel the error. After confirming that nothing is blocking the signal between the control unit and the foot control, wait for a while to check if the foot control is connected automatically. If it is not connected automatically, perform the following procedures.</p> <p>Turn the control unit power off and on again to connect to the foot control.</p> <p>Refer to "3.2 Installing foot control batteries" and replace the foot control batteries with new ones.</p> <p>Refer to "7.1.1 Pairing with foot control" and perform pairing.</p> <p>If this does not resolve the issue, please contact an authorized NSK dealer.</p>
"COM ERR LK" is displayed.	Link Module not connected	<p>Press the Cancel key to cancel the error. After confirming that nothing is blocking the signal between the control unit and the Link Module, wait for a while to check if the Link Module is connected automatically. If it is not connected automatically, perform the following procedures.</p> <ul style="list-style-type: none"> • Turn off and on both this product and VarioSurg 3 connected to this product. • Refer to the Link Module manual and try the pairing again. <p>If this does not resolve the issue, please contact an authorized NSK dealer.</p>
"COM ERR OS" is displayed.	Osseo 100+ not connected	<p>Press the Cancel key to cancel the error. After confirming that nothing is blocking the signal between the control unit and Osseo 100+, wait for a while to check if Osseo 100+ is connected automatically. If it is not connected automatically, perform the following procedures.</p> <ul style="list-style-type: none"> • Turn off and on both this product and Osseo 100+ connected to this product. • Charge the battery sufficiently with the Osseo 100+ charger and try the automatic connection again. • Refer to the Osseo 100+ manual and try the pairing again. <p>If this does not resolve the issue, please contact an authorized NSK dealer.</p>

8 Troubleshooting

Problem	Cause of problem	Action
"BATTERY FC" is displayed.	Foot control battery is low	Depress the speed control pedal of the foot control or press the Cancel key to cancel the error. Replace the batteries. If this does not resolve the issue, please contact an authorized NSK dealer.
"BATTERY OS" is displayed.	Osseo 100+ battery is low	Depress the speed control pedal of the foot control or press the Cancel key to cancel the error. Charge the battery sufficiently with the Osseo 100+ charger and try the automatic connection again. If this does not resolve the issue, please contact an authorized NSK dealer.

9 Spare Parts and Option List

Model	REF	Remarks
AC Power Cord	U260414	230V
AC Power Cord	U198013	230V: Only for Oceania
AC Power Cord	U371550	120V
Coolant Solution Hanger Post	U370152	
FC-81	Z1401001	Foot Control with Hanger
Battery Cover	Z1401068	For Foot Control.
SGL80M	Y1004211	Optic motor (with motor cord)
SG80M (2.0m)	Y1004212	Non-optic motor (with motor cord)
Protection Plug	10001595	
O-Ring	D0312457102	For protection plug.
O-Ring	D0312074080	For micromotor.
Tube Holder	E1198105	For micromotor.
Tube Holder	Y900083	Pack of 7
Calibration Bur	Z1057101	
Handpiece Stand	Z1402110	
Nozzle Holder	20000396	For X-SG20L.
Nozzle Holder	20000357	For SG20.
Internal Spray Nozzle	10000324	23mm
Cleaning Wire	20000512	
E-Type Spray Nozzle	Z019090	
Y-Connector	C823752	Used for branching the internal and external coolant irrigation.
Irrigation Tube	Y900113	Pack of 5 for 2.0m motor.
X-SG20L	C1003	Optic surgical handpiece.
SG20	C1010	Non-optic surgical handpiece.
X-DSG20L	C1068	Optic dismantling surgical handpiece.
X-DSG20	C1067	Non-optic dismantling surgical handpiece.
X-DSG20Lh	C1076	Optic dismantling surgical handpiece with hexagon.
X-DSG20h	C1075	Non-optic dismantling surgical handpiece with hexagon.
Link Module	Z1402001	For linking this product and VarioSurg 3.
Link Stand2	Z1452001	For installation of this product and VarioSurg 3.
VarioSurg 3	Y1002726	With foot control
VarioSurg 3 Non FT	Y1002248	Without foot control
Osseo 100+	Y1004176	Implant stability measurement system
Carrying Case	Y1004219	For components and accessories of this product. Only for Europe.
Carrying Case	Y1500395	For components and accessories of this product.
SG-CASE	S900040	Sterilization cassette. Only for Europe.

9 Spare Parts and Option List

Model	REF	Remarks
iCart Duo	S9090	For installation of this product and VarioSurg 3 and accessories. Only for Europe.
iCart L	Z299	For installation of this product and VarioSurg 3 and accessories.
PANA SPRAY Plus	Z182100	For high- and low- speed handpieces.

10 Specification

10.1 System specifications

■ Control unit

Model	NE335
Power Supply Voltage	AC100 - 240V
Frequency	50/60Hz
Power Consumption	Max. 240VA
Max. Pump Output	75ml/min
Dimensions	W245 x D235 x H90 mm
Weight	2.1 kg
Torque	5-80 N · cm



■ Micromotor

Model	SGL80M	SG80M
Rotation Speed Range	200 – 40,000 min ⁻¹ ±10%	
Input Voltage	DC36V	
Dimensions	Ø23.5 x L82.6 mm (without motor cord)	
Cord Length	2 m	
Optic	High CRI LED	-
Weight	0.2kg (with motor cord)	

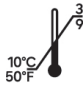■ Foot control

Model	FC-81
Electric Power Source	Three AAA/Micro/LR03 1.5V alkaline batteries
Dimensions	W260 x D185 x H65 mm
Weight	1.1 kg (with hanger)

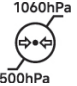■ Use, transportation and storage environment

	Temperature	Humidity	Atmospheric pressure
Use Environment			
Transportation and Storage Environment			

* No moisture condensation in the control unit.

* Use outside of these limits may cause malfunction.

Medical instruments that can be combined with the product:

- NSK VarioSurg 3
- NSK Osseo 100+

Refer to <Other products> in “2.2 Applicable products“.

10 Specification

- ISO 3964 (EN ISO 3964) compliant NSK Surgical handpiece

However, torque precision is not guaranteed for handpieces other than those listed in <Contra angle handpiece> in “2.2 Applicable products”.

The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by NAKANISHI INC. is under license. Other trademarks and trade names are those of their respective owners.

10.2 Bluetooth specifications

Frequency band	2.4GHz ISM band (2.402-2.480GHz)
Transmitting power	Class3 1mW
Modulation	GFSK
Channels	40channel 2MHz spacing
Compatibility	EN 300 328, EN 300 489-1, EN 301 489-17, EN 62479:2010

10.3 Classification of equipment

- Type of protection against electric shock:
 - Class I equipment
- Degree of protection against electric shock:
 - Type B applied part ⚡ (Applied part: handpiece)
- Method of maintenance (cleaning) and sterilization recommended by the manufacturer:
 - Refer to “5 Post-use Maintenance”
- Degree of protection against ingress of water as detailed in IEC 60529 (EN 60529):
 - Foot control: IPX8 (Protected against the effects of continuous immersion in water)
- Degree of safety of application in the presence of a flammable anesthetic mixture with air or with oxygen or nitrous oxide:
 - Equipment not suitable for use in the presence of a flammable anesthetic mixture with air or with oxygen or nitrous oxide.
- Mode of operation:
 - Non-continuous operation

10.4 Operation principle

Power is supplied to the control unit by operations on the foot control.
This makes the motor run and the bur attached to the handpiece rotate.

10.5 Factory settings (implant system)

The following are the initial factory settings. These are only reference setting conditions, so change them to the conditions for the implant system you use.

Implant system number	Program number		Maximum rotation speed [min ⁻¹]		Torque upper limit [N · cm]		Rotation direction [F/R]		Gear ratio [X:X]	Coolant solution flow level [0-5]		Light level [H/L]
1	1	5	2000	1200	-	-	F	F	20:1	3	3	H
	2	6	2000	25	-	45	F	F		3	3	
	3	7	2000	25	-	50	F	R		3	3	
	4	8	1600	25	-	20	F	F		3	0	
2	1	5	2000	800	-	-	F	F	20:1	3	3	H
	2	6	800	25	-	45	F	F		3	3	
	3	7	800	25	-	50	F	R		3	3	
	4	8	800	25	-	20	F	F		3	0	
3	1	5	800	400	-	-	F	F	20:1	3	3	H
	2	6	600	15	-	35	F	F		3	3	
	3	7	500	15	-	40	F	R		3	3	
	4	8	400	15	-	35	F	F		3	0	
4	1	5	800	300	-	-	F	F	20:1	3	3	H
	2	6	600	15	-	35	F	F		3	3	
	3	7	500	15	-	40	F	R		3	3	
	4	8	400	15	-	35	F	F		3	0	
5	1	5	1400	15	-	35	F	F	20:1	3	3	H
	2	6	1400	15	-	35	F	F		3	3	
	3	7	1400	15	-	40	F	R		3	3	
	4	8	800	15	-	35	F	F		3	0	
6	1	5	1400	1400	-	-	F	F	20:1	3	3	H
	2	6	1400	1200	-	-	F	F		3	3	
	3	7	1400	1000	-	-	F	F		3	3	
	4	8	1400	25	-	35	F	F		3	0	
7	1	5	1400	1000	-	-	F	F	20:1	3	3	H
	2	6	1400	15	-	35	F	F		3	3	
	3	7	1400	15	-	40	F	R		3	3	
	4	8	1200	15	-	35	F	F		3	0	
8	1	5	800	600	-	-	F	F	20:1	3	3	H
	2	6	600	15	-	35	F	F		3	3	
	3	7	600	15	-	40	F	R		3	3	
	4	8	600	15	-	35	F	F		3	0	

* Example) Implant system number: 3, Program number: 5

Maximum rotation speed: 400, Torque upper limit: -, Rotation direction: F,

Gear ratio: 20:1, Coolant solution flow level: 3, Light level: H

11 EMC Information (Electromagnetic Compatibility Information)

11 EMC Information (Electromagnetic Compatibility Information)

Guidance and manufacturer's declaration - Electromagnetic Emissions.

The product is intended for use in the electromagnetic environment specified below. The customer or the user of the product should assure that is used in such an environment. Guidance and manufacturer's declaration electromagnetic emissions.

Emissions test	Compliance	Electromagnetic environment - guidance
Radiated emissions CISPR 11/EN 55011	Group 1 Class B	The product uses RF energy only for its internal functions. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
Conducted emissions CISPR 11/EN 55011	Group 1 Class B	The product is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2/EN 61000-3-2	Class A (220-240V)	
Voltage fluctuations/flicker emissions IEC 61000-3-3/EN 61000-3-3	Complies (220-240V)	

Guidance and manufacturer's declaration - Electromagnetic Immunity

The product is intended for use in the electromagnetic environment specified below.

The customer or the user of the product should assure that it is used in such an environment.

Immunity test	IEC/EN 60601 test level	Compliance level
Electrostatic discharge (ESD) IEC 61000-4-2/EN 61000-4-2	± 8kV contact ± 2kV, ± 4kV, ± 8kV, ± 15kV air	± 8kV contact ± 2kV, ± 4kV, ± 8kV, ± 15kV air
Electrical fast transient/burst IEC 61000-4-4/EN 61000-4-4	Input a.c. power PORT ± 2kV Signal input/output parts PORT ± 1kV	Input a.c. power PORT ± 2kV Signal input/output parts PORT ± 1kV
Surge IEC 61000-4-5/EN 61000-4-5	Line to line ± 0.5kV, ± 1kV Line to ground ± 0.5kV, ± 1kV, ± 2kV	Line to line ± 0.5kV, ± 1kV Line to ground ± 0.5kV, ± 1kV, ± 2kV
Voltage dips IEC 61000-4-11/EN 61000-4-11	0% U _T ; 0.5cycle (0°, 45°, 90°, 135°, 180°, 225°, 270° and 315°) 0% U _T ; 1cycle and 70% U _T ; 25 (50Hz)/30 (60Hz)cycles Single phase ; 0°	0% U _T ; 0.5cycle (0°, 45°, 90°, 135°, 180°, 225°, 270° and 315°) 0% U _T ; 1cycle and 70% U _T ; 25 (50Hz)/30 (60Hz)cycles Single phase ; 0°
Voltage interruptions IEC 61000-4-11/EN 61000-4-11	0% U _T ; 250 (50Hz) /300 (60Hz) cycle	0% U _T ; 250 (50Hz) /300 (60Hz) cycle

11 EMC Information (Electromagnetic Compatibility Information)

Power frequency (50/60Hz) magnetic field IEC 61000-4-8/EN 61000-4-8	30A/m	30A/m
Proximity magnetic field IEC 61000-4-39	134.2kHz (Pulse modulation) 2.1kHz, 50% Duty, 13.56MHz (Pulse modulation) 50kHz 50% Duty	134.2kHz (Pulse modulation) 2.1kHz, 50% Duty, 13.56MHz (Pulse modulation) 50kHz 50% Duty

NOTE Ut is the a.c. mains voltage prior to application of the test level.

Guidance and manufacturer's declaration - Electromagnetic Immunity

The product is intended for use in the electromagnetic environment specified below.

The customer or the user of the product should assure that it is used in such an environment.

Immunity test	IEC/EN 60601 test level	Compliance level	Electromagnetic environment - guidance
Conducted RF IEC 61000-4-6/EN 61000-4-6	3 V rms 0.15 MHz - 80 MHz 6 V rms ISM bands between 0.15 MHz and 80 MHz	3 V rms 6 V rms	<p>WARNING: Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the product, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result. Interference may occur in the vicinity of equipment marked with the following symbol:</p> 
Radiated RF IEC 61000-4-3/EN 61000-4-3	3 V/m 80 MHz - 2.7 GHz	3 V/m	
	9 V/m in ISM bands 710 MHz, 745 MHz, 780 MHz, 5.24 GHz, 5.5 GHz, 5.785 GHz	9 V/m	
	27 V/m in ISM band 385 MHz 28 V/m in ISM bands 450 MHz, 810 MHz, 870 MHz, 930 MHz, 1.72 GHz, 1.845 GHz, 1.97 GHz, 2.45 GHz	27 V/m 28 V/m	

NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

11 EMC Information (Electromagnetic Compatibility Information)

- a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the product series is used exceeds the applicable RF compliance level stated above, the product should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the product series.
- b Over the 150 kHz to 80 MHz frequency range, the field strength should be less than 3V/m.

Cables and accessories	Maximum length	Complies with
Micromotor with Motor Cord	2.0 m	RF emissions, CISPR 11: Class B/ Group 1
AC Power Cord	2.0 m	Harmonic emissions (220-240V): IEC 61000-3-2 / EN 61000-3-2
Foot Control	non wired	Voltage fluctuations/flicker emission (220-240V): IEC 61000-3-3 / EN 61000-3-3
Link Module	non wired	Electrostatic discharge (ESD): IEC 61000-4-2/EN 61000-4-2
Osseo 100+	non wired	Electrical fast transient/burst: IEC 61000-4-4/EN 61000-4-4
		Surge: IEC 61000-4-5/EN 61000-4-5
		Voltage dips, short interruptions and voltage variations on power supply input lines: IEC 61000-4-11/EN 61000-4-11
		Power frequency (50/60 Hz) magnetic field: IEC 61000-4-8/EN 61000-4-8
		Conducted RF: IEC 61000-4-6/EN 61000-4-6
		Radiated RF: IEC 61000-4-3/EN 61000-4-3

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NSK Europe GmbH